

stryker[®]

Osteosynthesis

HA Apex[™] Pins

Pin Fixation System

**HA Coated Half Pins for Long
Term Fixation**




HA Apex™ Pins

Pin Fixation System

The HA coated Self Drilling/Self Tapping and Blunt Apex™ Pin range with different diameters, lengths and thread lengths provide the optimal choice to meet the needs of the treatment.


The HA coated threads reinforce the pin/bone interface.
HA Apex™ Pins are delivered sterile packed for direct use in the O.R.

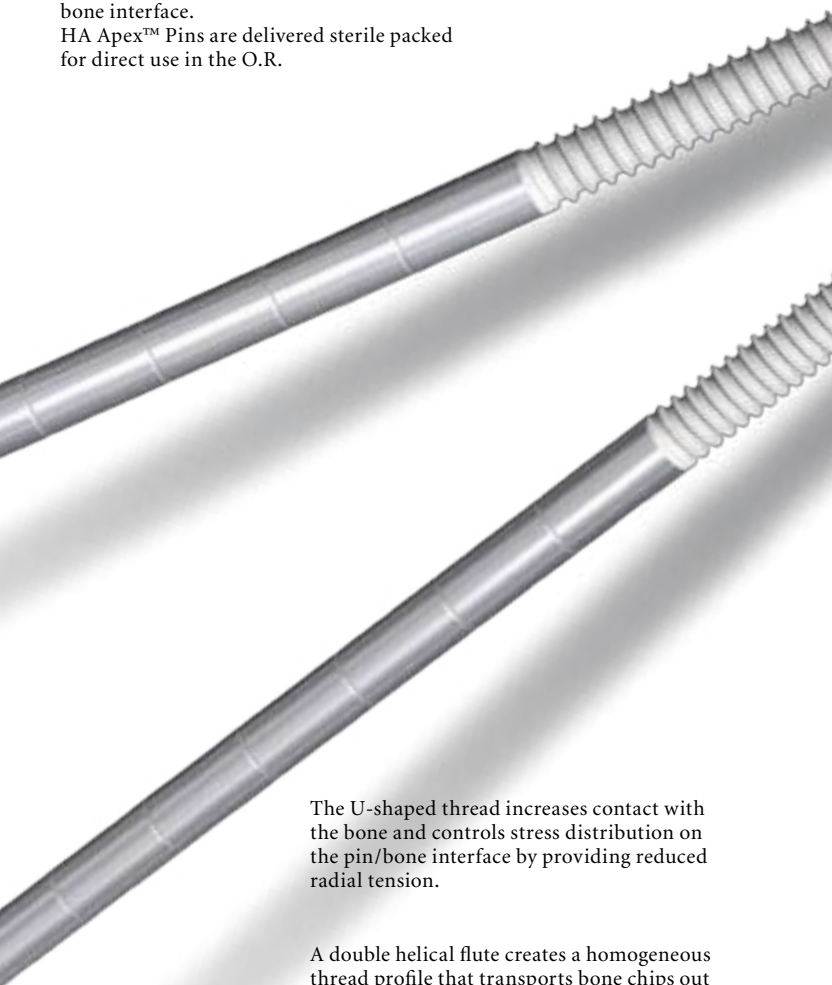

The U-shaped thread increases contact with the bone and controls stress distribution on the pin/bone interface by providing reduced radial tension.

A double helical flute creates a homogeneous thread profile that transports bone chips out of the drill hole for additional improvement of pin/bone interface.

The smooth tip design of the Blunt HA coated Apex™ Pin helps to reduce soft tissue irritation at the far cortex. Pre-drilling is required for this design.

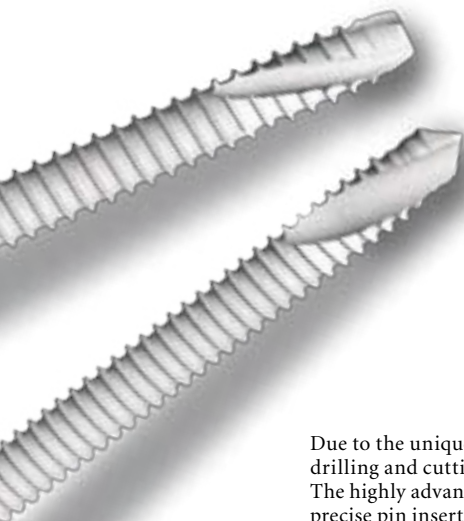

Due to the unique coating technology excellent drilling and cutting abilities are achieved. The highly advanced cutting geometry facilitates precise pin insertion with reduced insertion time and temperature, ensuring enhanced performance.

The cylindrical thread design improves bone purchase, and pull out resistance, and offers the possibility to back out the pin without compromising fixation.



Joint Replacements

Trauma, Extremities & Deformities

Craniomaxillofacial

Spine

Biologics

Surgical Products

Neuro & ENT

Interventional Pain

Navigation

Endoscopy

Communications

Imaging

Patient Handling Equipment

EMS Equipment

Stryker Trauma AG
Bohnackerweg 1
CH-2545 Selzach - Switzerland
www.osteosynthesis.stryker.com

The information presented in this brochure is intended to demonstrate a Stryker product. Always refer to the package insert, product label and/or user instructions before using any Stryker product. Surgeons must always rely on their own clinical judgment when deciding which products and techniques to use with their patients. Products may not be available in all markets. Product availability is subject to the regulatory or medical practices that govern individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area.

Stryker Corporation or its subsidiary owns the registered trademark: Stryker
Stryker Corporation or its subsidiary owns, uses or has applied for the following trademarks: Apex

Literature Number: 5075-4-002
LOT A2706

Copyright © 2006 Stryker
Printed in Germany

