

VariAx Elbow

Locked Plating System



Comprehensive Locking/Non-Locking Plating System

6 different anatomical plate designs of multiple lengths incorporate VariAx Polyaxial Locking Technology, allowing screws to be locked anywhere within a 30° cone.

180° & two Different 90° Distal Humerus Plating Options


Medial, Extended Medial, Lateral, Posterior Medial and Posterior Lateral plates allow for treating a broad range of elbow indications.

Anatomic and Contourable Design

Anatomically pre-contoured plates may be further adjusted to intra-operatively adapt to individual patient anatomy without damaging the locking mechanism.

Strength Where it is Needed

Plates are specifically designed with a bridging support between the diaphyseal and metaphyseal regions, making the plate stronger where the primary fracture line normally lies.

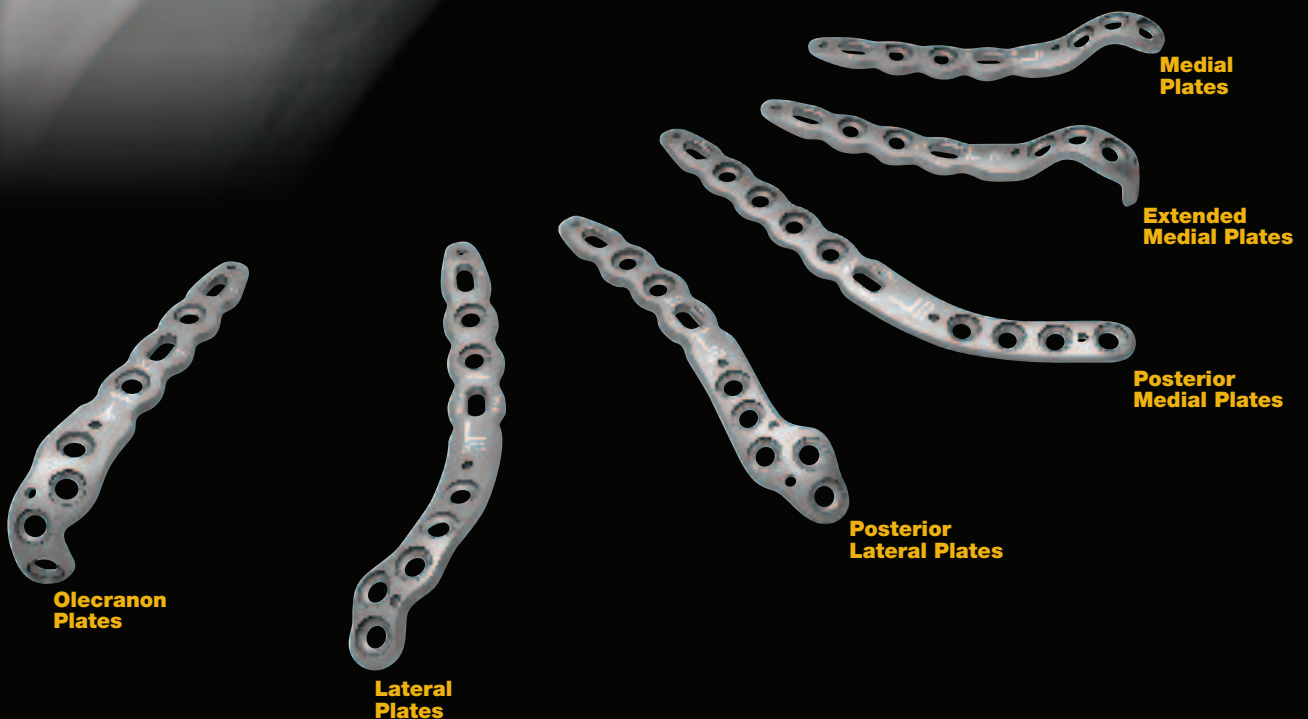

VariAx


Locked Plating Systems

VariAx is a comprehensive polyaxial locked plating system that provides the ability to treat various fracture patterns and deformities in a wide range of anatomical indications.


VariAx Distal Radius Locked Plating System




VariAx Hand Locked Plating System

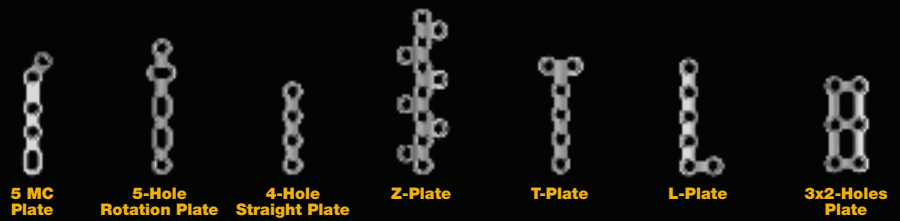

VariAx Foot Locked Plating System

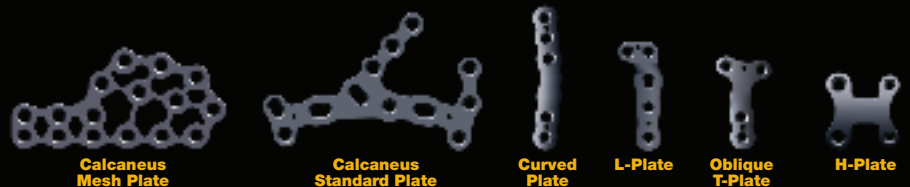

VariAx Fibula Locked Plating System

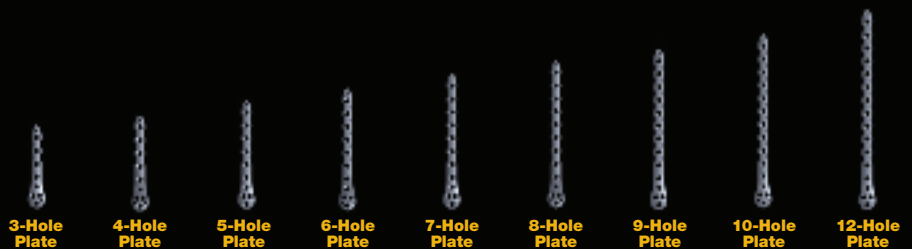

www.stryker.com

A surgeon must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that surgeons be trained in the use of any particular product before using it in surgery.

The information presented is intended to demonstrate the breadth of Stryker product offerings. A surgeon must always refer to the package insert, product label and/or instructions for use before using any Stryker product. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: Stryker, VariAx. All other trademarks are trademarks of their respective owners or holders. The products listed above are CE marked.

Literature Number: 982349
LOTA0910

Copyright © 2010 Stryker